

MN Sons Lab Solution

Laboratory Furniture and Fume Hood Manufacturer



Laboratory Furniture



Laboratory Fume Hood



1241

1. 11

11

Laboratory
Stand Alone Units

CONTACT US



 \bigcirc

+91 9825119622 / 7600824188





303/8, GIDC Industrial Estate, Makarpura, Vadodara. 390010.



ABOUT US

We are an ISO Certified laboratory furniture and fume hood manufacturer formerly known as Scientific Lab.

With our full in-house manufacturing capacity we bring you the same superior quality products with added services.

CERTIFCATIONS



• Laboratory Furniture Layout Planning.

- Product designed following the international standards.
- Customized Product for your unique requirement.
- Faster Delivery
- One year warranty on product and after sales service.



LABORATORY FURNITURE



Island Bench Pedestal Design



Types of Laboratory Furniture

- Wall Bench
- Island Bench
- Peninsular Bench
- Double Analytical Table

Wall Bench C-Frame Design

Accessories Available

Electrical Sockets
 6/16 Amp / Industrial Socket / BCH Metal
 Clad Sockets (Schneider or Legrand
 Make)
 Handles
 PP Black Handle / Aluminum Flush Handle
 / SS Handle(D Type)

• Hinges Ebco / Hettich

• Drawer Glides Ebco / Hettich

• Sink PP Sink / SS Sink / Ceramic Sink

• Table top Granite / Trespa / Metal

• Water Tap Three way / One way

• Peg Board PP / SS

• Cabinet Locks Godrej / Ebco / Hafele



Wall Bench Skirting Design

MANUFACTURING PROCESS

- Procurement of GI steel and CRCA steel from Tata Steels / Essar Steels.
- Sheet metal is cut into the needed component size on laser cutting machines.
- Bending of the components are done on semi automatic bending machines.
- The components are then spot welded together.
- Seven tank pretreatment process is performed to ensure that components are ready for powder coating.
- Powder coating is done by **electrostatic process** in spraying booth.
- Powder coated components are cured in the oven to obtain smoother finish. (70-80 microns)
- Components are assembled in semi knocked down condition and packed ready for dispatch.

LABORATORY FUME HOOD

Types of Laboratory Fume Hood

- Bench Mounted Fume Hood
- Low Bench Fume Hood
- Walk-in Fume Hood
- Customized Fume Chamber

The laboratory fume hood is the most common local exhaust ventilation system used in laboratories and is the primary method used to control inhalation exposures to hazardous substances. When used properly, fume hoods offer a significant degree of protection for the user. The Fume hoods are connected to the blowers mounted on terrace and through the PP Duct all the harmful gases are drawn out of the fume hood and released out of the laboratory. If the gases are too harmful to directly release in environment than scrubber system is used which filters out the toxic gases and make it safe to release it in atmosphere.




LABORATORY STAND ALONE UNITS



• Apron Cabinet



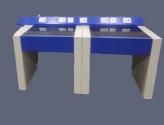





SS Stool



• Chemical Storage

• Lab Chair



• Anti Vibration Table

• SS Trolly



• Lockers

Cross over bench

OUR CLIENTS

